

**HYDRAULIC  
UNIVERSAL &  
HEAD-END CONTROLLED  
OPERATING TABLE**

**Model No.**

**ST-06**



**STURDY INDUSTRIAL CO., LTD.**

*MORE SAFETY - RELIABILITY - SATISFACTION*

**STURDY**  
APEX GROUP

# MULTI-FUNCTION OPERATING TABLE ST-06


## General Design

- \* Stable base construction with four castors for easy movement
- \* Base brake lock/unlock controlled by foot pedal
- \* Lateral tilt and trendelenburg adjusted by hand crank
- \* Base cover, side rails, accessories & clamps made of #304 stainless steel
- \* Table top elevation controlled by hydraulic system

- \* Table top was divided into 5 sections including head plate (U-shaped head rest optional), back plate, seat plate and leg plates
- \* 360° table top rotation
- \* With side rails for the insertion of X-ray cassettes (must equipped with acrylic sheet plates for X-ray radiography)

## Kidney surgery position

Table top can be adjusted to the Flex/Reflex positions for kidney surgery; controlled by head-end crank

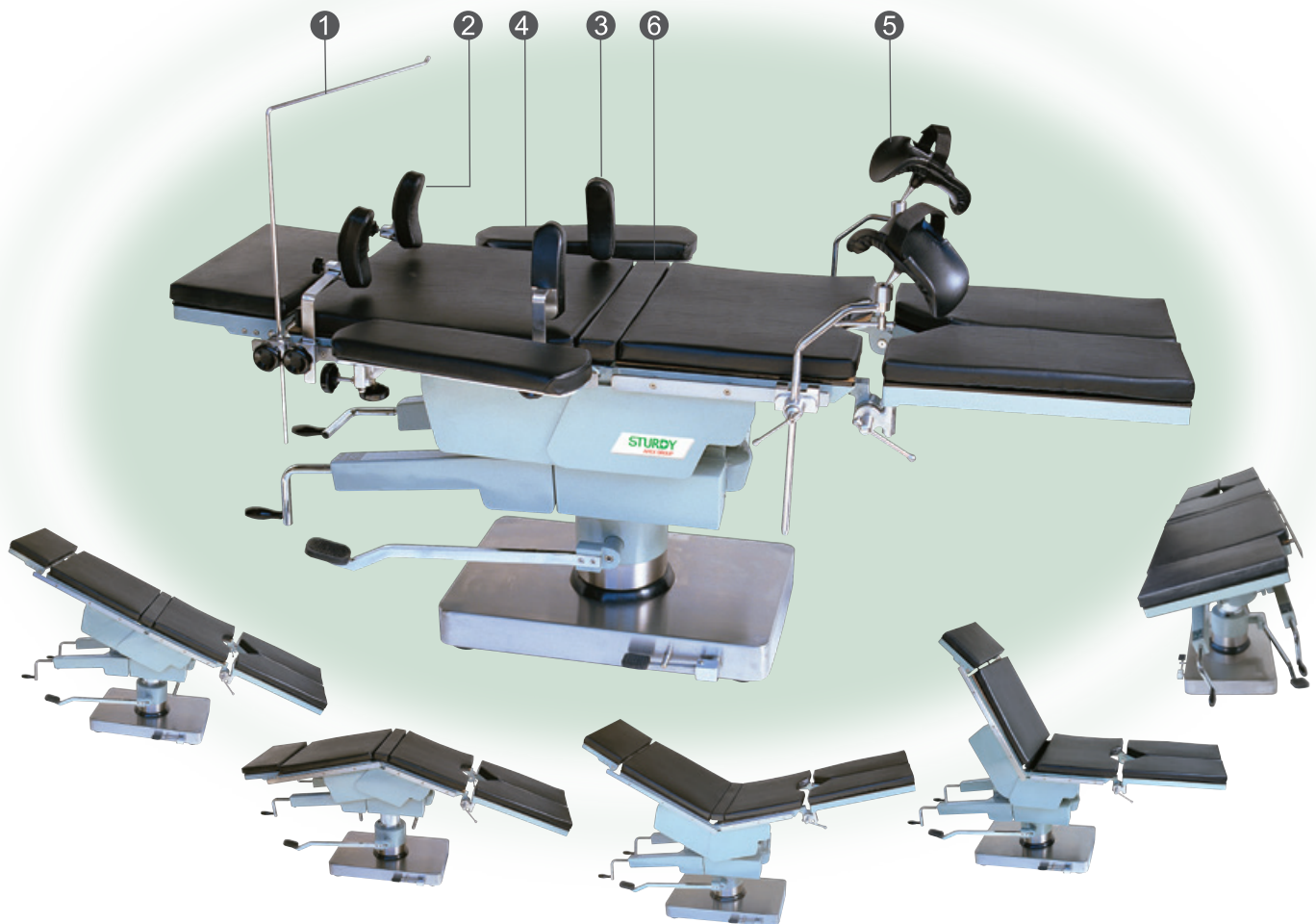

### SPECIFICATION:

# LENGTH OF TABLE TOP	200 CM
# WIDTH OF TABLE TOP	51 CM
# TABLE TOP LOW / HIGH	68 - 98 CM
# LATERAL TILT LEFT / RIGHT	25° / 25°
# TRENDLENBURG / REVERSE	22° / 28°
# BACK SECTION UP/DOWN	80° / 20°
# SEAT SECTION UP/DOWN	20° / 20°
# HEAD SECTION UP/DOWN	30° / 90°
# LEG SECTION UP/DOWN	20° / 90°
# BRAKE FUNCTION BY FOOT PEDAL	

### STANDARD ACCESSORIES:

1. SCREEN FRAME (TAC-01)	1 SET
2. SHOULDER SUPPORTER (TAC-10)	1 PR.
3. BODY SUPPORT (TAC-11)	1 PR.
4. ARM BOARD (TAC-14)	1 PR.
5. LEG CRUTCHES (TAC-13)	1 PR.
6. MATTRESS	1 SET

### OPTIONAL ACCESSORIES:

# RESTRAIN STRAP (TAC-02) / PC.
# HAND GRIPS (TAC-08) / PR.
# U-SHAPED HEADREST (TAC-03) / SET
# RAISED ARMREST (TAC-04) / SET
# X-RAY RADIOGRAPH TABLE TOP (TAC-12) / SET
# CONTAMINATION BUCKET (TAC-09) / SET

\* STURDY reserves the right to improve the equipment shown without prior notice.